

RANGE OF POSTURAL AIDS FOR BED

# SYST'AM® P902T - P904T - P906T / HEEL PAD

HEEL PAD MADE OF VISCOELASTIC FOAM WITH MULTIBEARING SURFACE AND MEMORY EFFECT

BED RELATED

## MATERIALS

80 kg/m<sup>3</sup>



Visco foam



Removable POLYMAILLE® integral cover



NF EN ISO 597 - 1 & 2

Foam maintenance :



(Do not immerse in water)

Cleaning cover :



## INDICATIONS

- Used as prevention, the heel pad to be positioned at the bottom of the bed is particularly well suited to cases of moderate loss of mobility of the lower limbs with a moderate to very high risk of sores forming.
- When used to assist in treatment, this pad enables treatment of the sore in cases of low to normal mobility.

### SYST'AM® P902T / HEEL PAD (without cover)



### P904T / HALF-HEEL PAD

Enables different positions or pressure relief for each lower limb.

### FABRIC COATED WITH BI-STRETCH POLYURETHANE


Reduces friction and shear effects.  
Supple and soft to the touch (comfortable).  
Favours the exchange of gases (steam, sweat) :  
- fights against maceration.  
Impermeable material:  
- better hygiene,  
- longer support system lifespan.  
Washable at 90°C, can be decontaminated using cold sprays.  
Treated to resist fire.  
With a non-slip lower face to help to stay in place.  
In multi-patient use, it is preferable to buy one new cover per patient.

### SYST'AM® P902T / HEEL PAD (with cover)



### NON-SLIP LOWER FACE

Helps the system to stay in place.

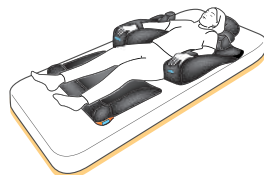

DESIGNATION	ITEM CODE	SIZE L x W x H (cm / inches)
P902T / Heel Pad	P902T1HW	73 x 64 x 10,5 cm / 28,7 x 25,2 x 4"
P904T / Half-Heel Pad	P904T1HW	73 x 32 x 10,5 cm / 28,7 x 12,6 x 4"
P906T / Heel Pad with abduction wedge	P906T1HW	73 x 64 x 25,5 cm / 28,7 x 25,2 x 10"



WARRANTY



03/2016

