



# SYST'AM® P141M / POLYMULTI - in 3 sections

A MOULDED WAFFLE MATTRESS MADE OF HR FOAM WITH HIGH DENSITY

## MATERIALS

40 kg/m<sup>3</sup>



HR foam



PUR POLYFILM fitted bed sheet



POLYMAILLE® integral cover



### Foam maintenance:



(Do not immerse in water)

### Cleaning POLYFILM cover:



### Cleaning POLYMAILLE® cover:



## INDICATIONS

### Risk




P141M / POLYMULTI Mattress can be used, when:

- 1- Confinement in bed is necessary resulting from diseases or handicap and there is at the same time an increased risk of developing pressure sore according evaluation scales (Braden, Norton, Waterlow, ...).
- 2- Pressure sore (Stages 1 & 2 according EPUAP) has developed.



17 cm

→ A sole-type base giving greater stability.

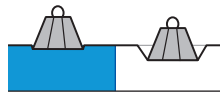

### PERFORMANCE OF THE MULTIBEARING CONCEPT


→ The multibearing surface enables a better pressure distribution through the transfer of supports from areas with low bearing capacity toward areas with high bearing capacity:

#### WHITE AREAS:

Areas with low bearing capacity (softer foam) granting an increased prevention on areas at risk (sacrum, heels, shoulder blades).

#### BLUE AREAS:


Areas with high bearing capacity (firmer foam) for parts of the body showing a lower level of risk.

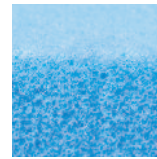

### PERFORMANCE AND COMFORT OF THE SYST'AM® HIGH RESILIENCE FOAM

- Because it is moulded and of high density (40 kg/m<sup>3</sup>), it ensures durable flexibility, comfort and a soft touch.
- It offers as well a better resistance against external wear and tear, for a longer lifespan of the mattress.

### POLYMULTI®: Mattress with multibearing surface

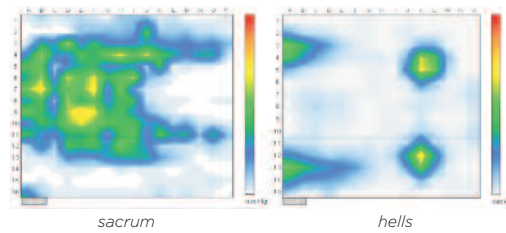

sacrum

heels

## AVAILABLE COVERS






POLYFILM FITTED BED SHEET FOR SINGLE PATIENT USE

- A breathable and hyper-extensible hypoallergenic polyurethane film.
- Unfailing adaptability: Unlike PVC, this material is extremely strong and offers lasting resistance against wear and tear.
- POLYFILM is totally impermeable to liquids.



POLYMAILLE® INTEGRAL COVER FOR MULTI-PATIENT USE

- Hypoallergenic polyurethane enduction (anti-bacterial, anti-germ, anti-fungal).
- POLYMAILLE® is bi-extensible, to reduce tear and shear forces.
- Zip over 3 sides. Non-slip base.

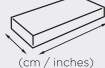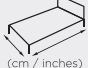DESIGNATION	ITEM CODE	 (cm / inches)	 (cm / inches)	 min - MAX (kg / lbs)
P141M / POLYMULTI with POLYFILM fitted bed sheet	P141MMPFDHNC	86 x 197 / 34 x 78	90 x 200 / 35 x 79	35 - 130 / 77 - 286
P141M / POLYMULTI with POLYMAILLE® integral cover	P141MMPMHNC	86 x 197 / 34 x 78	90 x 200 / 35 x 79	35 - 130 / 77 - 286

AVAILABLE ONLY NOT COMPRESSED



WARRANTY 2 YEARS