



# SYST'AM® P914L / UNIVERSAL CUSHION

POSTURAL AID AND BEDSORE PREVENTION

## MATERIALS

80  
kg/m<sup>3</sup>



Visco foam



Removable  
POLYMAILLE®  
integral cover



NF EN ISO 597 - 1 & 2

### Foam maintenance:



(Do not immerse in water)

### Cleaning POLYMAILLE® cover:



## INDICATIONS


- Universal positioning help.



### VISCOELASTIC FOAM WITH SHAPE MEMORY


- Defined body moulding and increased body surface area in contact with the foam.
- More effective distribution of pressure: reduction of transcutaneous pressure on highrisk zones, facilitated blood flow.
- Increased comfort and more effective maintenance of patient position.

### AVAILABLE VERSIONS

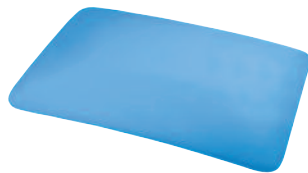

P914L / UNIVERSAL CUSHION

- This is the perfect supporting aid: it helps support the patient, relieves the pressure on the heels, protects the internal lateral knee areas, ankles...
- This cushion can also be used as postural help in Rest chair, wheelchair and shell seat.
- Facilitates positioning for care staff.
- Used by itself or with other devices.

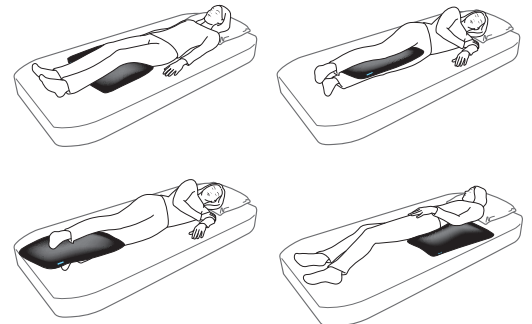

## FEATURES OF THE COVERS

### FABRIC COATED WITH BI-STRETCH POLYURETHANE

- Reduces friction and shear effects.
- Supple and soft to the touch (comfortable).
- Favours the exchange of gases (steam, sweat):
  - fights against maceration.
- Impermeable material:
  - better hygiene,
  - longer support system lifespan.
- Washable at 90°C, can be decontaminated using cold sprays.
- Treated to resist fire.
- In multi-patient use, it is preferable to buy one new cover per patient.



DESIGNATION	ITEM CODE	SIZE L x W x H (cm / inches)
P914L / Universal cushion	P914L1HW	60 x 40 x 13 cm / 24 x 15,7 x 5"

